

PMR-30 platform rocker – fixed tilt

Compact but highly functional fixed-angle platform rocker in the Grant bio range, providing a smooth side to side rocking motion for gentle sample agitation in tubes, culture flasks, dishes and boxes.


- **Variable speed: 5 to 30 oscil/min**
- **Fixed 7° tilt angle**
- **Load up to 0.5 kg**
- **Continuous or timed operation, with automatic switch-off**



Reliable and extremely quiet motor produces regulated and reproducible rocking throughout the speed range

Compact, sturdy construction with a low profile and small footprint – fits neatly into the smallest workspace

Simple timer setting, with large easy-to-read LED display indicating clearly the time remaining

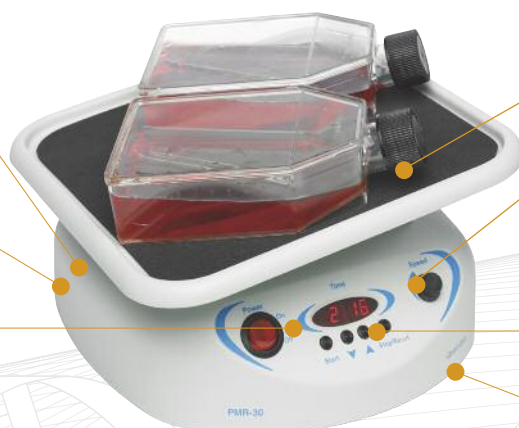

Smooth, non-slip mat supplied as standard – prevents vessels from slipping

Simple graduated speed setting – from 5 to 30 oscil/min

Select either continuous or timed operation – the integral electronic timer ensures accurate count-down for repeatability of time-sensitive incubations


Low voltage cord easily fits through incubator door gaskets

Applications: ■ Life science - ideal for minigel destaining after electrophoresis, conducting reactions of Northern, Southern and Western blots, immunoblots, hybridisation washes, cell culture, visualisation of blots, immunostaining, protein electrophoresis.

PMR-100 platform rocker – adjustable tilt

Adjustable angle and speed platform rocker, providing soft or intensive side to side rocking for optimal mixing of samples. Suitable for use in cold rooms or incubators with ambient operating temperature range +4°C to 40°C.

- **Variable speed: 1 to 99 oscil/min**
- **Tilt angle range up to 10°**
- **Load up to 5 kg**
- **Continuous or timed operation, with automatic switch-off**

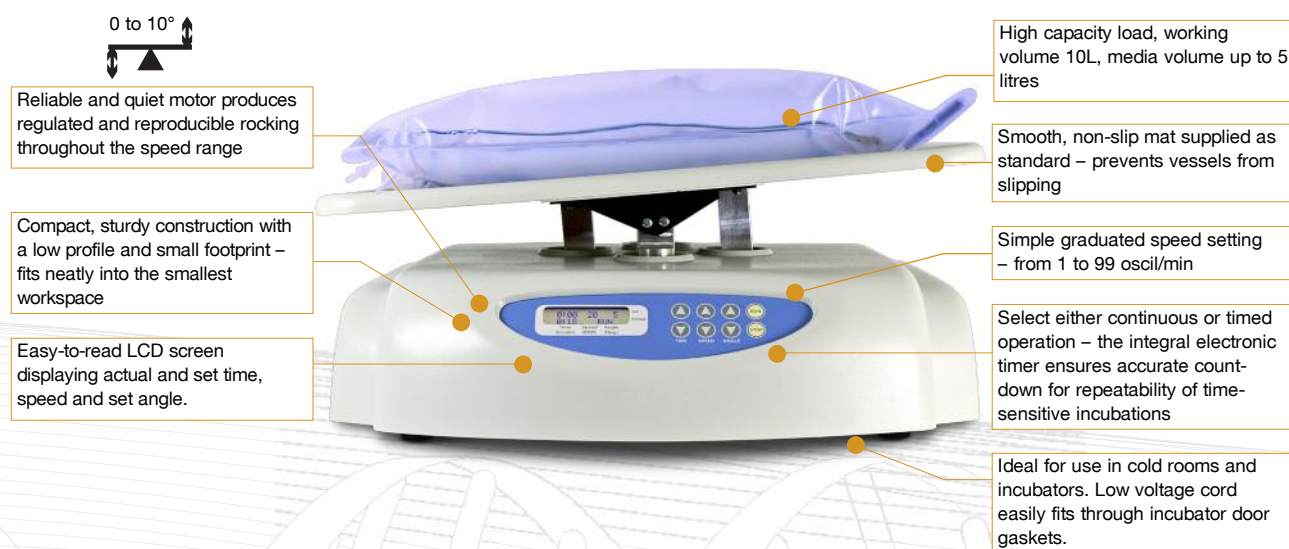



Applications: ■ Life science - ideal for minigel destaining after electrophoresis, conducting reactions of Northern, Southern and Western blots, immunoblots, hybridisation washes, cell culture, visualisation of blots, immunostaining, protein electrophoresis.

Rockers and rotators » specifications

Grant-bio compact rockers and rotators – models and specifications

● = optional

		PMR-30	PMR-100
		Fixed tilt platform rocker	Adjustable tilt platform rocker
		 h: 120 mm d: 205 mm w: 220 mm weight: 2kg	 h: 250 mm d: 480 mm w: 400 mm weight: 13kg
Speed	oscill/min	5 to 30	1 to 100
Fixed tilt angle	°	7	-
Adjustable tilt angle	°	-	10 (from 51-99 oscill/min)
Tilt angle range	°	-	0-9 (increment 1°) (from 1-50 oscill/min)
Timer, with automatic switch-off		1 min to 23hrs 59mins	1min to 99hrs 59min
Maximum continuous operation time		168 hrs	
Platform dimensions (working area)	mm	210 x 210	460 x 360
Maximum load	kg	0.5	5
Display		4 digit LED	LCD, 16 x 2 character
Ambient temperature range	°C	4 to 40	4 to 40
Input voltage	V dc	12	12
Input current	A	0.5	0.5

Accessories

PDM – dimpled mat



●

-



Thank you for reading this data sheet.

For pricing or for further information, please contact us at our UK Office, using the details below.



UK Office

Keison Products,

P.O. Box 2124, Chelmsford, Essex, CM1 3UP, England.

Tel: +44 (0)330 088 0560

Fax: +44 (0)1245 808399

Email: sales@keison.co.uk

Please note - Product designs and specifications are subject to change without notice. The user is responsible for determining the suitability of this product.